Deta
Elis
Holding

Gesundheit
Innovationen
Mensch

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Deta
Elis
Holding
Health
Innovations
Human

Device for enhancing your own
defensive systems
and correction of natural
electromagnetic radiation of
natural biosystems

DeVita
Energy

Purpose of the device

The device “DeVita Energy” is a generator of weak electromagnetic radiation, identical to the radiation of natural organic systems.

Designed to enhance your own protective systems and correct natural electromagnetic radiation of natural biosystems. The principle of corrective action is based on the use of frequencies according to the methods of R. Rife, R. Voll, F. Morrel.

Attention! The device may only be used by children under adult supervision. Children are prohibited from playing with the device. Persons with physical disabilities and intellectual disabilities may only use the device under the supervision of appropriate personnel.

Box set

“DeVita Energy” device.....	1 pcs.
Passport/User’s Manual.....	1 pcs.
Packaging.....	1 pcs.
USB charging cable.....	1 pcs.
Neck strap.....	1 pcs.

Technical characteristics

1. Number of programs.....	11
2. Frequency range, Hz.....	0,1-15000
3. Accuracy error of frequency setting no more than.....	$\pm 0,003$
4. Time of continuous work no less than, hours.....	15
5. Keyboard type.....	sensor
6. Screen type, pixels.....	TFT 128x160
7. Power supply:	
- Li-Po battery 3,7V, built-in, pcs.....	1
8. Maximum charging time, hours.....	6
9. Overall dimensions not more than, mm.....	103x45x16
10. Device’s weight no more than, kg.....	0,1
11. Average service life not less, years.....	5

Safety information

Device is:

- Safe to use
- Non-toxic
- Not flammable
- Not explosive
- Not a source of dangerous radiation

When charging the Class III battery, devices must be powered by an Exclusively Safe Low Voltage AC adapter.

Environment protection

The built-in battery contains environmental pollutants. The battery must be removed at the end of its useful life and before sending it for disposal. The battery must be disposed of at a designated collection point. 

Instructions for care, storage and transport

Device without packaging must be stored in a room temperature from 10 to 35 °C and a relative humidity of 80%.

The device must be protected from direct moisture affection and active chemicals.

Do not disassemble the device.

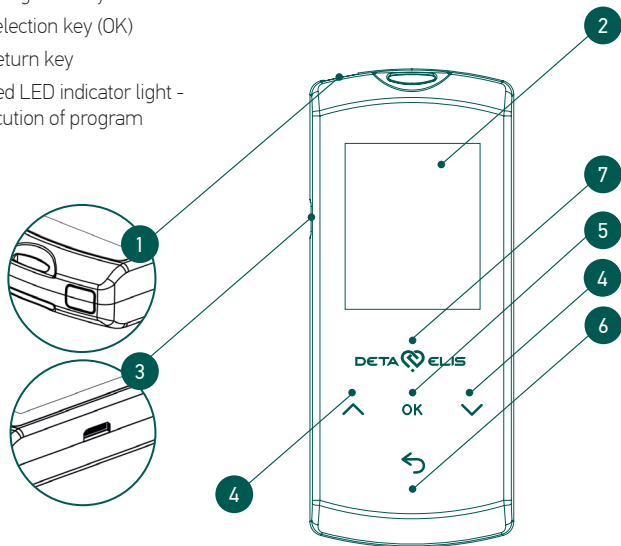
Do not drop the device.

Do not drop other objects onto the device.

We recommend keeping the packaging during the warranty period.

Device's visual description

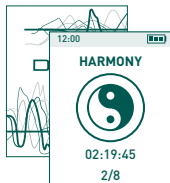
1. ON/OFF button
2. Display
3. Charging port
4. Navigation keys
5. Selection key (OK)
6. Return key
7. Red LED indicator light - execution of program

Pic. 1

Working with the device

To turn the device on press ON/OFF button **1**. Screen timeout is 15 seconds from the last press on any button. Display may be turned off before that time by pressing ON/OFF button **1** again. If the display went out before you entered the required commands, it must be turned on by pressing the ON/OFF button **1** again.

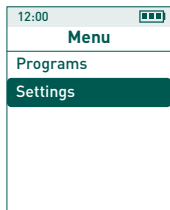
After pressing the **1** ON/OFF button (Pic.1) you will see the home page, and after 2 sec. the device automatically enters the program selection mode.

To switch to the main menu mode, press the return key.

Working with the main menu

Submenu "Settings"

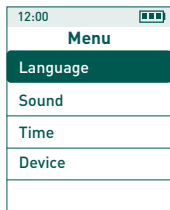
Choose "Settings" using the navigation keys

Press the selection key

Submenu "Settings" will be shown on the display

Language selection

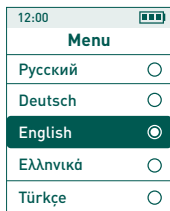
Choose "Language" using the navigation keys

Press the selection key

Submenu "Language" will be shown on the display

Choose the preferred language with the navigation keys

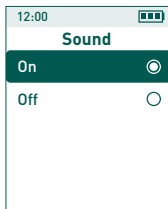
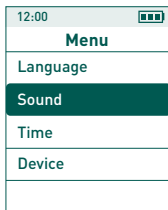
Press selection key

To switch to the device main menu mode, press the return key

Attention! The default language is English. If there is no need for changing the language you can skip this step.

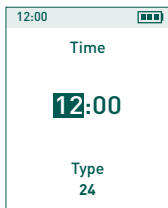
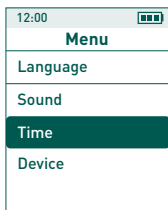
Sound settings

-  Choose "Sound" using the navigation buttons
-  Press the selection key

Submenu "Sound" will be shown on the display

-  Choose between "On" or "Off" with the navigation keys
-  Press the selection key
-  To return to submenu "Settings" press the return keys

If you select "Off", the sound of pressing any key is turned off.

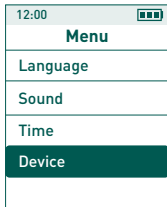
Time settings

-  Choose "Time" with the navigation keys
-  Press the selection key

Submenu "Time" will be shown on the display

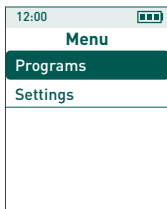
-  Use the navigation keys to set hours, minutes and type
-  Use the selection key to select "Hours", "Minutes" or "Type"
-  To return to submenu "Settings" press the return key

Device

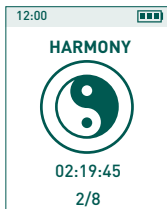
-  Choose "Device" with the navigation keys
-  Press the selection key

Manufacturer serial number and program version will be shown on the display.

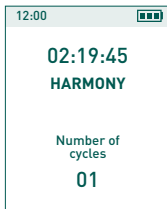
Submenu "Programs"

-  Choose "Programs" with the navigation keys
-  Press the selection key

Program selection

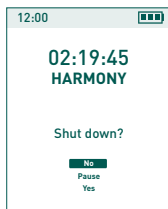
-  Choose the necessary program (eg Harmony) with the navigation keys
-  Press the selection key

Cycles set screen will be shown on the display

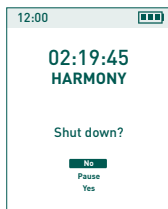
-  Set the needed number of cycles with the navigation keys
-  Press the selection key

The program will start, which will be signaled by a flashing red LED indicator light pic 1 (7).

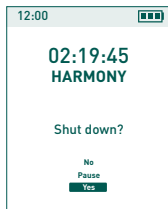
Submenu with the ongoing program will be shown on the display

Timer in the upper part of the screen shows the countdown till the program finishes.

Information about the program in progress

In case you will need the information about the program in progress press **1** ON/OFF button (Pic.1). Display will show the submenu of the program in progress.

Interruption/cancellation of the program

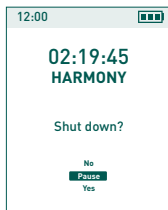
If the program needs to be cancelled

Choose "Yes" with the navigation keys

Press the selection key

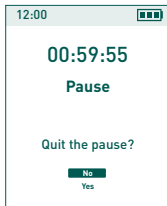
If the program needs to be interrupted (paused)

Choose "Pause" with the navigation keys

Press the selection key

Display will show the "Pause" submenu
If you want to quit the pause

Choose "Yes" with the navigation keys

Press the selection key

If you will not choose "Yes" in "Quit the pause?" submenu, your device will be turned off automatically within an hour.

Attention! When the device is paused, the red LED indicator light fig. 1 (7) will flash at a slower rate!

Using the device with the launched program

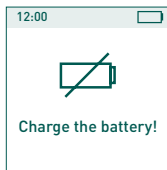
After starting the program, place the device with the it's back facing your body, for example, in your pocket.

Charge the device

Battery charge indicator

If the battery indicator is white-colored and two  or three  internal lines are visible – than your battery is charged.

If the battery indicator is red-colored, one internal line is visible  – your battery is running low.

In case of low charge or full discharge of the battery, connect the micro-USB connector of the cable to the charging port of the device, and the USB connector to the USB port of the working computer. The battery indicator will begin to signal the charging process.

How to charge the battery

In case of low charge or full discharge of the battery, connect the micro-USB connector of the cable to the charging port of the device, and the USB connector to the USB port of the working computer, or an external charger with output 5V-1A. The battery indicator will begin to signal the charging process.

Manufacturer's warranty

The manufacturer guarantees full compliance with the technical specifications of the device "DeVita Energy", provided that the user strictly observes the instructions for operation, transport and storage.

The warranty period for the device is 24 months from the date of sale through the retail network.

In the absence of the date of sale and the stamp of the trade organization in the passport, the warranty period is calculated from the date of production of the device by the manufacturer.

During the warranty period, the owner has the right for free repairs. Warranty repairs are carried out at service centers authorized by the manufacturer. Delivery of the defective device to the place of repair is carried out at the expense of the buyer.

If there is a violation of external covers and / or in the absence of the certificate warranty, no claims are accepted and no warranty repair is performed.

The following faults are not covered by the warranty:

- defects due to improper use;
- defects caused by natural disasters;
- in case of violation of the integrity of seals;
- in case of presence of external defects (cracks, fractures etc.)

The buyer has the right to replace the defective device with a new one during the warranty period in the following cases:

- the device is repaired for the fourth time
- the device is unrepairable

Attention! The manufacturer reserves the right, without prior notice, to make changes in the design of the device and software without affecting the existing characteristics.

Methodical recommendations

For the first time, the idea of the presence of energy in the human body and the dependence of health on its circulation was expressed by the ancient healers of China, India, Egypt. Today it has been scientifically proven that every human organ emits certain frequencies. Under the influence of unfavorable external factors, full of stresses of modern life, poor ecology, its frequencies deviate from their natural parameters, as a result of which is depression, apathy, and chronic insomnia occur. The device „DeVita Energy“ emits harmonious electromagnetic frequencies that equalize the bioenergetic indicators of a person, trigger the body's regulation functions based on resonance phenomena, and reduce bio- and geopathogenic loads. „DeVita Energy“ gently and safely affects the human body, without disturbing its energy balance either as a whole or locally. This method of influence on the body is based on the developments of Franz Morel, the author of the method of bioresonance therapy. Correctly created 8 programs of electromagnetic influence restore and strengthen physiological indicators, relieve psycho-emotional stress, energize, give a good mood and excellent health, harmonize a person at all levels of his being.

Programs of „DeVita Energy“

„Abundance“ program. Working time – 2: 17: 14.

Description: people always have certain limiting beliefs, fears that do not make it possible to enjoy the financial resource or properly dispose of it. In this regard, fears arise: the fear of losing and the fear of having money (they can also be stolen, and if they are easily obtained, then a certain guilt complex arises). All this can cause intense tension, anxiety and other emotional experiences that are reflected in our body. The program helps to align the attitude towards money, regardless of a person's financial capabilities. In addition, it affects the change in attitudes towards material goods and partnerships, promotes the disclosure of creative abilities.

Application: if necessary, the program can be used daily, combining with the programs: „Relationship“, „Health“, „Harmony“.

„Harmony“ program. Working time – 02: 19: 45

Description: the program has a beneficial effect on the nervous system, balances the psycho-emotional state, helps to remove hidden blockages, subconscious negative beliefs, fears. It helps to improve sleep, unleash full potential, which

allows you to achieve goals in the best way, sets you up for positive.

Application: the second major wellness program, used daily in conjunction with the „Health“ program.

„Protection“ program. Working time – 01: 42: 01.

Description: working with people you can often experience negative feelings: headache, aggression, depression, loss of appetite, discord in partnerships, etc. All this can lead to functional disorders.

The program creates an excellent barrier against negative energy and does not allow it to harm health and mood. It will be useful for whose specialties are associated with work in large teams and in the stream (doctors, trainers, teachers, psychologists, social workers).

Application: the program is used as needed.

„Health“ program. Working time – 02: 59: 23.

Description: the program is aimed at maintaining the body at a functional level, restoring the energy potential of damaged organs, and rejuvenating.

Application: the first major wellness program, used daily in conjunction with the „Harmony“ program.

„Slimness“ program. Working time – 01: 01: 36.

Description: the program is aimed at eliminating psychosomatic causes that are factors in weight gain. Our subconscious mind can provoke overeating and, as a result, the deposition of extra pounds in response to psycho-emotional stress, which is caused by conflicts in the family, at work, household and emotional dissatisfaction. In addition, the program will be useful for people who want to gain weight, since it gives a feeling of inner calmness, quiet joy, which contributes to the process of proper food intake and its complete assimilation.

Application: if necessary, the program can be used daily in combination with the „DeVita Energy“ programs: „Relationship“, „Health“, „Harmony“.

„Relationship“ program. Working time – 01: 16: 11.

Description: used to harmonize relationships - in the family, between partners, friends, to resolve social problems. This program expands the possibilities of communication, mutual understanding, tolerance and love, helps to eliminate conflicts. Gives a sense of unity with the world. For introverts, it will be useful in order to release internal "clamps" and unleash creativity.

Application: if necessary, the program can be used daily in any combination with other „DeVita Energy“ programs.

„Longevity-5G“ program. Working time – 01: 15: 01.

Description: the program helps to maintain the general level of health at good values, contributes to the protection of DNA from the negative effects of electromagnetic radiation from household appliances, the Internet, etc. Adjusts the body according to its reference „sound“ (solfeggio frequency). As a result, the amount of vital energy, clarity of mind, consciousness increases, a mood for a healthy lifestyle appears, and creative abilities are activated.

Application: if necessary, the program can be used daily in combination with programs: „Health“, „Harmony“.

Attention! In the first days of the program, there may be: slight dizziness, nausea, drowsiness.

„Energy“ program. Working time – 00: 41: 12.

Description: the program energizes, gives strength for further development, has a healing effect, protects against bio- and geopathogenic factors.

Application: if necessary, the program can be used daily in any combination with other „DeVita Energy“ programs.

„Energy Pro“ program. Working time – 00: 09: 19.

Description: The program is aimed at restoring lost energy, increasing immunity and body tone. It helps to feel a burst of strength, increase the vital potential necessary for the realization of your desires.

Application: The program can be used daily in the morning and afternoon as needed.

„Anti age“ program. Working time – 03: 51: 05.

Description: To slow down aging, preserve youth and energy is a cherished dream of many people in different periods of history, from antiquity to the present day. To maintain a youthful face and body, women and men have resorted to various methods. „Anti age“ program is another, but the most modern method of rejuvenation, which starts the process of tissue regeneration, removal of toxins from cells, recovery of the body after 50 years of age. The program has an impact not only on the physical, but also on the psychological state, allowing you to feel young, active and confident for many years.

Application: The program can be used daily in combination with Health, Harmony, „Longevity-5G“ programs.

„Inner Peace“ program. Working time – 00: 25: 12.

Description: Due to various factors, actions or events of the present time, people are increasingly experiencing emotional tension, stress, accumulation of negative emotions, which negatively affects the state of health. The program helps to relieve mental tension, recover from stress, support emotional and physical well-being, give a sense of inner peace and harmony with the world.

Application: The program is used if necessary. In case of strong emotional tension, stress, it is allowed up to 5-7 repetitions in a row or alternation with the program “Harmony” until the improvement of well-being.

Recommendations for using the device „DeVita Energy“

The device is recommended for everyone who has to experience:

- The harmful effects of the city and the environment
- Physical and emotional fatigue
- Chronic fatigue syndrome
- Great physical activity
- Depression, stress, anxiety, irritability
- Dizziness, weakness, drowsiness

Effect after using the device „DeVita Energy“:

- Restores and enhances physiological performance
- Relieves psycho-emotional stress
- Removes psychosomatic blocks
- Charges with energy and positivity
- Launches the body's self-regulation functions
- Regulates drainage and detoxification functions of the body
- Reduces bio- and geopathogenic loads
- Protects against electromagnetic smog
- Rejuvenates

Conclusion

The Wellness concept – healthy eating, moderate exercise, fresh air, giving up bad habits – largely reduces the risk of various diseases. However, in the modern world it is not always possible to consistently adhere to this concept. This is especially true for megalopolises with their constant stress, fast-paced lifestyle and dysfunctional environment.

„DeVita Energy“ devices are designed to help the body restore lost energy, increase immunity and maintain vitality at a high level.

The higher the energy of a person, the better his health, the more potential and opportunities he has in order to create his life the way he wants.

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EC Declaration of Conformity Nr. W14

The company **DETA-ELIS EUROPA GMBH** with a seat at Kaiserstr. 74, 63065 Offenbach, Germany

Herewith to DECLARE

That the equipment described below complies with all applicable at the time of construction essential requirements of the **Directive 2011/65/EU (RoHS)** and the **Directive 2014/30/EU (EMC)**.

This declaration relates exclusively to the equipment in the state in which it was placed on the market, and excludes alterations or manipulations, components which are added and/or operations carried out subsequently by the final user without our official consent.

Equipment: DeVita
Type / Model: Energy
Standards Applied: EN 55014-1:2006+A1:2009+A2:2011
EN 55014-2:1997+A1:2001+A2:2008
EN 61000-3-2:2006+A1:2009+A2:2009
EN 61000-3-3:2008

Date: 23 January 2018

Authorized Representative: The General Manager Ms. Konopleva Tatyana

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CONCLUSION: IN REFER TO SGS TEST REPORT NO.: H0FA0005
NO DEVIATION FROM THE STANDARDS HAVE
BEEN IDENTIFIED. THEREFORE THE TESTED
PRODUCT IS IN COMPLIANCE WITH THE
APPROPRIATED STANDARDS.

This Verification of Compliance has been granted to the applicant based on the results of tests, performed by Laboratory of SGS-CSTC Standards Technical Services Co., Ltd. on sample of the above-mentioned product in accordance with the provisions of the relevant specific standards. The CE marking as shown below can be affixed, under the responsibility of the manufacturer, after completion of an EC Declaration of Conformity and compliance with all relevant EC Directives.

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SGS Germany GmbH, Hofmannstr. 50, D-81379 Munich is testing facility for
ELECTRO-MECH. COMPONENTS TESTING ENVIRONMENTAL SIMULATION ELECTROMAGNETIC COMPATIBILITY
MEDICAL DEVICE TESTING BATTERY TESTING PRODUCT SAFETY TELECOM CONFORMANCE TESTS

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Member of the SGS Group (Société Générale de Surveillance)

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Deta-Elis Europa GmbH